



Overview:

A very strong and versatile electric height adjustable podiatry treatment chair with dual leg rests with extending heel pads. The leg sections manually lower from +10 to -90° for various lower limb procedures. The backrest has manual adjustment from +85° to -15° and the fold-down armrests can be folded back for patient transfer or even removed if required. The chair can also be used for other lower limb treatment procedures.




Benefits

- Lifetime frame & 5 year parts guarantee - see website for further details
- Electric operation of chair height from 52 to 101cm
- Adjustable angle backrest for supporting the patient in sitting position
- The backrest is fitted with fold-down armrests - removable if required
- Positive and negative adjustment of the backrest angle from -15° to +85°
- The foot section angle is fully adjustable from +10° to -90°
- The heel pads can be extended to position the foot for various procedures
- Safe working load and lifting capacity of 250kgs (550lb / 39 stone)
- Retractable wheel system provides easy positioning of the chair
- Low patient transfer height of 52cm (20")
- Wide choice of options, colours and accessories to meet user requirements

Overall Dimensions

Length: 176cm (69") approx
 Width: 63cm (25") approx
 Height: 52cm - 101cm approx
 (20" - 40") approx


Upholstery Dimensions

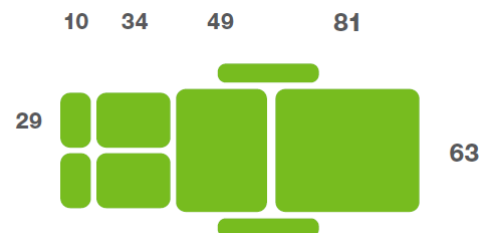

Options & Accessories:



Matching Operators Chair

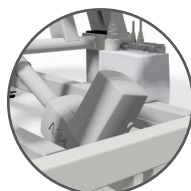

Battery Back-up



Motorised Backrest



Foot Switch Operation



D Cushion




Paper Roll Holder

NHS Approved Supplier



t: +44 (0)1449 767887
 e: sales@plinthmedical.com
 www.plinthmedical.com

Operational Parameters

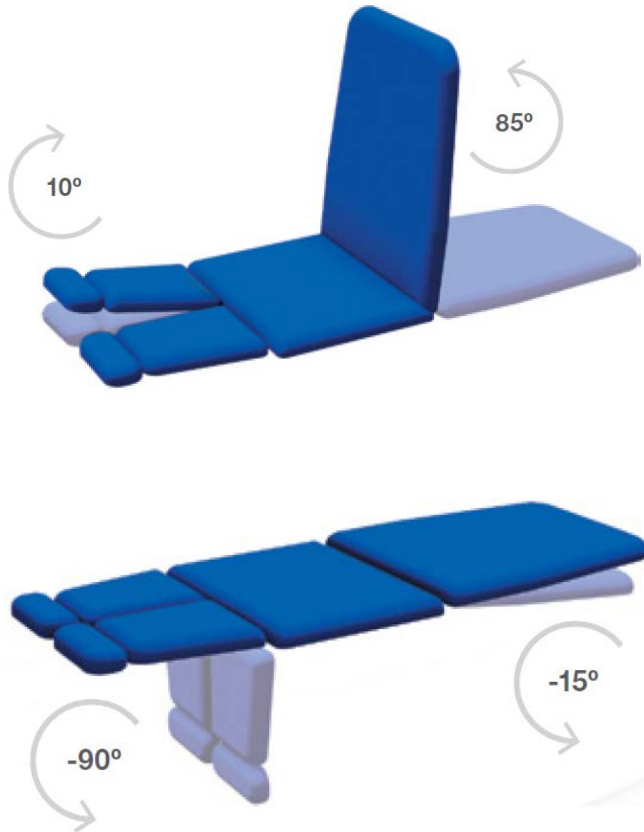

Model 93CD - Electric Version

Standard Features

- Designed and manufactured in the UK
- Lifetime frame & 5 year parts guarantee
- 16 year average lifespan
- Electric height operation
- Comfort foam
- Easy operation
- Low running costs
- Hardwearing anti-microbial vinyl
- Safe working load 250kg
- Retractable castors
- Hoist clearance 130mm
- Divided & extending legs
- Removable fold back arm rests

Technical Standards

- IP Rating - IPx4
- EEC Classification - MDD93/42/EEC Annex IX Class 1
- Electrical Standard - IEC601 Part 1 + 60601 P (EN)
- Vinyl Flame Retardant to CRIB 5 (BS 7176)
- Foam to BS 5852: Part 2: 1982, Ignition Source 5 [CRIB 5]
- CE Marked
- Externally audited to ISO9001:2015 by LRQA

Delivery Promise

- Delivery to UK mainland in just 7-10 working days
- Faster service may be available if required
- Overseas delivery 7-10 days plus shipping time

Upholstery Vinyl

- Flame retardant to BS 5852 Ignition Source 5 - CRIB 5 (BS 7176) Medium Hazard
- Independently tested to verify effectiveness at reducing the spread of MRSA, E.coli & other similar bacteria & fungi
- Rated at 500,000 Martindale rubs without damaging or marking the top layer, making it suitable for all high-traffic, heavy-duty areas
- Phthalate free
- Please note: If the upholstery colour is not specified, 'Almond' will be supplied



Plinth Medical Ltd
Barric Lane
Occold, Eye
Suffolk, United Kingdom
IP23 7PX

t: +44 (0)1449767887
e: sales@plinthmedical.com
www.plinthmedical.com

Plinth | Medical